

# **Prescription bariatric wheelchair**

Designed for users up to 39 Stone (250kg) 39 stone (250kg)

The Condor is an aluminium wheelchair with a unique folding mechanism which allows it to carry 39 Stone (250kg).

> The Condor is a highly adjustable wheelchair that can be configured to suit the user.

> > Adjustments include; seat to ground, rear wheel position, backrest angle, armrest height and seat depth.

> > > It is also a great looking wheelchair why not try it at your local Karma retailer?



Folding frame






 Flip back and detachable armrests

Prescription based seat and back system.

- Seat Widths 22", 24", 26", 28", 30" Back heights A: 16", 17", 18" B: 18", 19", 20"
- Back angles adjustment -5°, 0°, +5°, +10°, +15°
- 2 Seat Depths 18", 20", 22", 24"
- Seat to Ground Heights 15.5" ~ 19.5"
- Adjustable rear wheel position Rear wheels height and depth adjustable
- Adjustable footplates

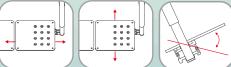

• Footrests are a multi-directional adjustable footplate.

KARMA Mobility Limited





#### Specifications

Model	Width	F/R wheels		<u>Gr</u>		- +	<u><u></u></u>	<u>Gr</u>	<u>F</u>		<u>Br</u>		Max. User Weight
KM-BT10	22"	200×50/24"	56	46/51	78	46	25	118	47	46/48.5/51	97/99.5/102	30.2	250
	24"	200×50/24"	61	46/51	83	46	25	118	47	46/48.5/51	97/99.5/102	31.0	250
	26″	200×50/24″	66	46/51	88	46	25	118	47	46/48.5/51	97/99.5/102	31.8	250
	28″	200×50/24"	71	56/61	93	46	25	123	49.5	46/48.5/51	99.5/102/104.5	33.5	250
	30″	200×50/24″	76	56/61	98	46	25	123	49.5	46/48.5/51	99.5/102/104.5	34.3	250

\* KARMA reserve the right to modify information herein without further notice.

±1cm ±0.5kg

## 🔥 Features

- 1.Maximum user weight 39 Stone (250kg).
- 2.Folding aluminium frame for transport and storage.
- 3.Adjustable rear wheel position.
- 4.Detachable flip back armrests.
- 5.Ergonomic armpads.
- 6.New footrest hanger release mechanism.
- 7. Adjustable footplates.
- 8.Slowing and parking brakes.
- 9.Bracing bar that can be used to push the wheelchair.
- 10.Unique anti tip / tipping lever.
- 11. Highly adjustable seat and backrest.

## 🔥 Options

#### **Options accessories:**

- Quick release rear wheel
- Height adjustable armrests
- Elevating legrests

#### **Miscellaneous Accessories:**

- IV Pole
- Oxygen Bottle Holder

Sold and serviced by:

104-405A3 Rev.00 OCT. 2011

### KARMA Mobility Limited

www.karmamobility.co.uk / e-mail: info@karmamobility.co.uk Jnit 6, Target Park, Shawbank Road, Redditch, B98 8YN U.K. FEL: 0845-630-3436 FAX: 0845-630-3736