



Agile

The Agile is simple, but never ordinary. It is built with aircraft grade aluminum, which makes it light in weight and durable.





Easy Transfer

The flip-back armrests make moving from side to side easy and quick.

The 80° removable footrests ensure a safe transfer.



Smooth and Agile Movement

Agile wheelchair is designed with the carer in mind.

With advanced engineering, this chair is very easy to push and manoeuvre.



Compact and Lightweight

The Agile weighs from 12.7 KG complete, which is lighter than your carry on luggage.

With the backrest and seat folded, the chair's compact size makes it very convenient to carry and travel with.

Lightweight and Transportable

The Agile is simple, but never ordinary. Agile is built with aircraft grade aluminum, which makes it lighter in weight (12.7KG) and durable in use.



The Curved Line Design helps guide the seat upholstery into a nice and neat fold. Quick-Release Rear Wheels allow lighter carrying weight and smaller storage space.



▲ Black Frame, 18" Seat Depth, 24" Rear Wheel With Folding Backrest

The Agile is available in a range of configurations, that will fit the body dimensions of most users.

Seat Height

505mm

Seat Width

16" 、 18" 、 20"

Seat Depth

410mm/ 460mm



Backrest Frame Types

Foldable

Rear Wheel Types

Fixed/ Quick-released


Quick to Adapt. Ready for Use



Attendant-Propelled Model



Self-Propelled Model

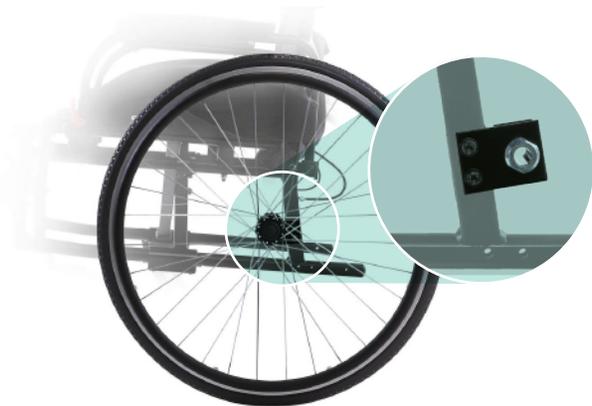

Adjustable Rear Wheel Axles

By adapting the optional rear-wheel extension plate, you can configure the position of rear wheel axles and optimise Agile's maneuverability or stability.

Robust and Reliable

The Agile is approved for transportation in accordance with ISO 7176-19 Crash Test. This means it is suitable to be used as a seat in a motor vehicle for transportation.





Optional Accessories

Tension Adjustable Backrest

Unlike traditional backrest designs, the tension adjustable backrest is made up of high-quality Velcro straps. By changing the tension of each strap, the ideal curvature can be created to fit the user's back.



Attendant Brakes



Flip-back,
Height Adjustable
Armrests



Elevating
Legrests



Anti-Tips



Seat Belt



Crutch Holder



I.V. Pole Holder



Oxygen Bottle
Holder



Frame Colour



Flat Black

Upholstery Colour



Black



Scan to learn more about the Agile.

Agile Specifications

Seat Width* Seat Depth (inch)	16*16	18*18	20*18
Seat Width (mm)	430	480	530
Seat Depth (mm)	505	555	555
F/R wheels (inch)	F: 8 R: 14F/ 24Q	F: 8 R: 14F/ 24Q	F: 8 R: 14F/ 24Q
Seat Height (mm)	505	505	505
Backrest Height (mm)	510	510	510
Armrest Height (mm)	220	220	220
Calf Length (mm)	380-480	380-480	380-480
Max. User Weight (kg)	125	125	125
Total Weight (kg)	12.7-14.2	13.0-14.4	13.4-14.9
Overall L x W x H (mm)	Min. 1025*600*1015 Max. 1065*605*1020	Min. 1085*650*1015 Max. 1120*655*1020	Min. 1085*650*1015 Max. 1120*655*1020
Overall Folded L x W x H (mm)	Min. 750*280*760 Max. 850*280*750	Min. 800*280*760 Max. 900*280*750	Min. 800*280*760 Max. 900*280*750

Remark : Karma reserves the right to modify information herein without further notice.

Unit:cm ±1cm ± 0.5kg

Rev.01 Mar.2021



KARMA MOBILITY LIMITED

TEL : 0845-630-3436

FAX : 0845-630-3736

Email : info@karmamobility.co.uk

www.karmamobility.co.uk

soma

Designed and Made by Karma