



UNLIMIT YOUR WORLD

We at Küschall believe anything is possible when you have the power to achieve it. The Küschall K-Series offers modularity and performance in a single product. Unlimit your Strength, Unlimit your Freedom, **Unlimit your world**.

KÜSCHALL K-SERIES

K-series is an **icon** of design, a trendsetter in the monotube frame construction. The ingenious **modular** seating system will fit your needs and its **dynamic** set up will offer you a superb driving performance, to keep you on the move.

HYDROFORMING

The hydroformed front frame will bring the design of your K-Series to another level. The new technology brings more stiffness and rigidity, for an unrivalled driving performance, adding striking lines and a clear design language to a timeless icon.

HYDROFORMED FRAME

 Ergonomic grip, ideal for repositioning and transfers, frame angle available in both 75° and 90°.

LIGHT ACTIVE BRAKE

The new aluminium brake, with a minimalist design, brings lightness to unprecedented levels. Precision even in the smallest details.

TAPERED FRAME

 Tapered frame allows the wheelchair to contour to your body for a sleek aesthetic.

FOOTBOW

 Completely redesigned with a flat surface for greater comfort, and perfect positioning.

Küschall® ĸ-series

NEW BACKBRACE

 Ergonomic back brace for increased grip when lifting during transfers. Lower brace positioning, ideal for Matrx backrests.

BACKREST

 Folding backrest with six possible angle adjustments for different positioning needs.

SEAT MODULE

At the heart of the K-series. The seat module offers five different centre of gravity positions, allowing a choice of high activity or increased stability within the wheelchair, as standard.

CASTOR HOUSING

New aluminium profile welded to the front frame tube. Consistent design throughout the entire wheelchair.

FORK

BBTER

New clean design and integrated bearings. Reduces the amount of debris collecting in the bearings, for increased performance.

HI IFT

DEUVLE

K

KEY FEATURES

CUSTOMISATION

ICONIC

With its iconic minimalist open frame and lightweight design, the K-Series gives you the ability to enjoy everyday activities; whether it's getting to and from work, or pursuing your active hobbies, the K-Series gives you the freedom to be you.



MODULAR

Adjustments come as standard on the K-Series, so finding your perfect position is always guaranteed!

Possible Adjustments:

Front seat height Rear seat height Centre of gravity Backrest height and angle

DYNAMIC

All adjustments on the K-Series are on the seat module. This will ensure your castor fork angle is always fixed at the correct angle, ensuring a sharp and responsive performance for a superior driving experience.



HYDROFORMED FRONT FRAME

Hydroforming technology has added even more personality to the K-Series. An Iconic look that today becomes even more recognisable, with its straight edge design upfront. Complemented by its increased strength and rigidity in the frame, the K-Series driving experience is more responsive than ever.

CHOOSE THE FRAME COLOUR




2 CHOOSE THE DECAL STYLE





3 CHOOSE THE DECAL COLOUR

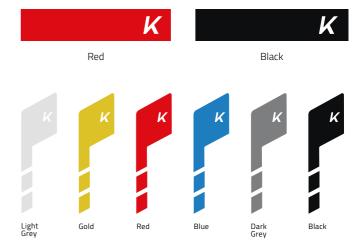

TECHNICAL SPECIFICATION

SEAT WIDTH 320 – 500 mm in increments of 20 mm



SEAT HEIGHT

front: 450 – 520 mm rear: 390 – 490 mm in increments of 10 mm



CAMBER

0°/3°/6°



SEAT DEPTH 350 – 525mm in increments of 25mm



FRAME ANGLE 75° / 90° V-frame or tapered



TOTAL WEIGHT UNLOADED approx. 7.8 kg



LOWER LEG

LENGTH

290 – 480 mm

in increments of 10 mm

5

TOTAL

WIDTH

SW + 170 mm

TRANSPORT WEIGHT approx. 5.1 kg



BACKREST HEIGHT 270 – 480 mm in increments of 15 mm



TOTAL LENGTH 75°: approx. 830 - 1055 mm 90°: approx. 760 - 985 mm



MAX. USER WEIGHT 130 kg



BACKREST ANGLE 74°/78°/82°/ 86°/90°/94°



CENTRE OF GRAVITY

80 mm / 95 mm / 110 mm / 125 mm / 140 mm



CRASHTEST

Approved - Fulfills the requirements of ISO 7176-19



All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration. Weight and measurements can change with the different configurations. For further information please check the user manual on our website.



UNLIMIT YOUR WORLD

Invacare Limited

Pencoed Technology Park - Pencoed Bridgend CF35 5AQ - United Kingdom Tel: +44 1656 776 222 Fax: +44 1656 776 220 E-mail: uk@invacare.com Sales Order E-mail: ordersuk@invacare.com www.invacare.co.uk Invacare Ireland Ltd Unit 5 - Seatown Business Campus Seatown Road - SWORDS - County Dublin - Ireland Tel. +353 1 8107084 Fax +353 1 8107085 Email: Ireland@invacare.com www.invacare.ie