Invacare® | Leo®





Safe, solid, secure: the Invacare Leo

The *Invacare Leo* is a scooter designed for individuals who value their independence and wish to get out and about unaided. Available in 3 and 4 wheel options, safety is a key feature of the *Leo*, but this does not detract from its stylish and sporty looks. *Leo* offers individuals the freedom and confidence to enjoy their essential daily outings and leisure excursions. With high standards of design and safety engineering, the *Leo* sets a new market standard.





For a safe, powerful drive

The *Invacare Leo* offers a smooth, safe drive in a scooter that is easy to control and manoeuvre. The *Leo's* intelligent lighting package includes indicators that automatically switch off and a rear brake light for safe driving. In addition to a visual warning for the battery level display, there is an audible warning (beep) when the battery needs recharging.

Compact

The **Leo** is a compact mini-scooter that can be easily dismantled for transportation during longer trips or weekend breaks. Easy manoeuvrability for getting around corners in the local shopping centre, museum or leisure complex becomes refreshingly easy. Access to the battery and electronics for maintenance or servicing is much easier with the **Leo** as only one part needs to be removed.


Reliable

Individuals can trust the **Leo** to take them wherever they need to go. From a trip to the shops or a ride to the park, the **Leo** is sturdy and reliable thanks to its exceptionally high quality electronic engineering and power. These elements are protected from corrosion and

water, (splash protection) which safeguard the vital workings of the scooter against dirt and water. A two-step disengaging lever protects the scooter from free-wheeling if inadvertently knocked out of drive mode. The lighting system is cleverly designed to keep functioning even when one light fails.

......

Features and Accessories

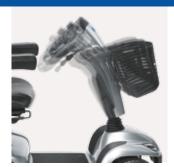

Finite tiller movement
Individuals can easily adjust
the tiller to suit their needs.



Brake lightRear brake lighting warns other people that the scooter is slowing down, thus

preventing collisions.



Comfort seatingSwivel functionality and sliding rails ensure optimal seating position.



Two-step disengaging leverThis unique concept prevents the scooter from accidentally freewheeling.

Features and Accessories





Built in anti-splash guardsProtects electronics and the transaxle from dirt and water.



Disassembly into five partsNo tools are required to dismantle *Invacare* **Leo**



Crutch & cane holderKeeps your crutch or cane stable during the drive.



Lockable front boxSecure box for personal valuables.



Seatbelt



Large mirrorCan be mounted left or right.




Heavy-duty storage cover Keeps the scooter clean and dry.



Charging off-boardSimply remove batteries to charge with the off-board charger kit.



Seat suspension moduleAdds extra comfort during driving.

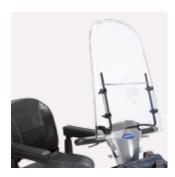

WindscreenProtects the individuals from the wind.

Technical data













4 wheels	470 mm	410 mm
3 wheels	470 mm	410 mm


460 - 560 mm

475 mm 475 mm

590 mm 590 mm

1200 mm



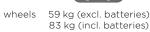



2 x 36 Ah







4 wheels 54 kg (excl. batteries) 3 wheels

136 kg 78 kg (incl. batteries)

2 x 36 Ah

2620 mm 2320mm

60 mm 70 mm









4 wheels 240 W 38 km 3 wheels 10° * 240 W 36 km 8 km/h

Colour of shrouds







Onyx Blue



......

Ruby Red -Only available in 4 wheel version

Upholstery



Black vinyl

You Tube



Invacare Limited Pencoed Technology Park -Pencoed - Bridgend CF35 5AQ -United Kingdom Tel: +44 1 656 776 222 Fax: +44 1656 776 220 E-mail: uk@invacare.com Sales Order E-mail: ordersuk@invacare.com www.invacare.co.uk

Invacare Ireland Ltd Unit 5 - Seatown Business Campus Seatown Road - SWORDS -County Dublin - Ireland Tel. +353 1 8107084 Fax +353 1 8107085 Email: Ireland@invacare.com www.invacare.ie

© 2016 Invacare International Sàrl All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Leo - UK - 07/2016 - 1503875





Follow us on YouTube:

Invacare Ltd - UK



^{*} Seat in maximum fifth rearward locking position from the front.