



Invacare® | **Matrix® Flo-tech™ Contour Lo-Back**

Allows the individual to sit back further into their seat

**matrix**  
SEATING SERIES

**INVACARE**  
Yes, you can.®

The *Invacare **Matrix Flo-tech Contour Lo-Back Cushion*** has been designed with a lower rear contour allowing the individual to sit further back into seat. The additional contouring and support under the thighs stabilises the individual further, thus reducing the damaging effects of shear.



### Contour Lo-Back Cushion optimises support surface contact area

Comparison with standard wheelchair cushion, showing improved posture through use of the **Matrix Flo-tech Contour Lo-Back Cushion**.

## Features and Options



### Ingress resistant cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for transportation.



### Available in 79 sizes

The Invacare **Matrix Flo-tech Contour Lo-Back Cushion** is available in a range of 79 sizes, from 14-24" widths and 13-22" depths. (See size chart for details).



### Contoured seating

Offers improved stability and postural support. Castellations reduce shear and friction forces on tissue.



### Ultra-fresh™ added to all cushions for antimicrobial and odour protection







Antimicrobial treatment has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria and helps provide longer lasting freshness.

### Moulded Foam Cushion

Moulding allows **Matrix Flo-Tech** cushions to provide a three-dimensional form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature	Warranty
					
A range of 79 sizes	103 mm	1.38 kg	152 kg (24 stone)	80°	3 years


**Matrix Flo-tech Contour Lo-Back Cushion**

Manufactured to European Standards:  
BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)


1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

## Colours

Black Smooth

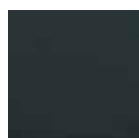

Black wipedown

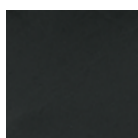

Black spacer fabric

