



*Invacare®* | ***Matrx® Flo-tech™ Contour***

Contoured moulded foam for improved posture

**matrx®**  
SEATING SERIES

**INVACARE**  
Yes, you can.®

The surface of the *Invacare* **Matrix Flo-tech** Contour cushion is contoured into sections. Its design supports the client and increases pressure redistribution. The combination of improved posture and increased support surface area will contribute significantly to pressure ulcer prevention. This deep cushion is suitable for heavier people.

## Moulded foam cushion

Moulding allows **Matrix Flo-tech** cushions to provide a three-dimensional form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity.

## Contoured seating surface

Offers improved stability and postural support; castellations reduce shear and friction forces on tissue.

## Features and Options



### Ingress resistant cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.



### Available in 84 sizes

The *Invacare* **Matrix Flo-tech Contour** is available in a range of 84 sizes, from 14-24" width and 14-22" depths. (See size chart for details).



### Integral sag rigidizer option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.



### Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria and helps provide longer lasting freshness.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

	Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature	Warranty	
<b>Matrix Flo-tech Contour</b>	A range of 84 sizes	103 mm	1.38 kg	152 kg (24 stone) 191kg (30 stone) on sizes 22" width/depth and above	80°	3 years	Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-1 Match, and BS 7176 (Crib 5)
							1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

## Colours


Black Smooth



Black wipedown



Black spacer fabric

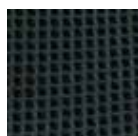

Baltic Blue

