



Invacare® | **Matrix® Vi**

Moulded foam with visco insert for protection and positioning

matrix®
SEATING SERIES

INVACARE
Yes, you can.®

The *Invacare **Matrx Vi Cushion*** is designed to provide skin protection and positioning for individuals at moderate to high risk of skin breakdown. The unique three part foam construction targets skin protection where it is needed most, and features excellent stability, support and comfort.



① **Outer cover**

Moisture resistant and breathable. Includes slip resistant base, hook and loop fasteners, lifting strap.

② **Inner cover**

Moisture resistant. Easily wipes clean. Zippered cover allows access for cushion modifications.

③ **Top layer**

Soft (ILD) High Resilient (HR) foam for optimal comfort and support. Pre-ischial contour provides stability to pelvis. Abductor/adductor contour promotes midline leg positioning.

Made with Ultra-Fresh™ for antibacterial and odour protection.


④ **Visco foam insert**

Moulded visco-elastic foam ischial relief insert. Anatomically targeted to ischial/sacral area. Provides immersion to sensitive bony prominences.


⑤ **Bottom layer**

Firm High Resilient (HR) foam for support and durability. Trochanteric shelf provides stability to hips. Waterfall front edge. Rear radius for sacral support.

Features and Options

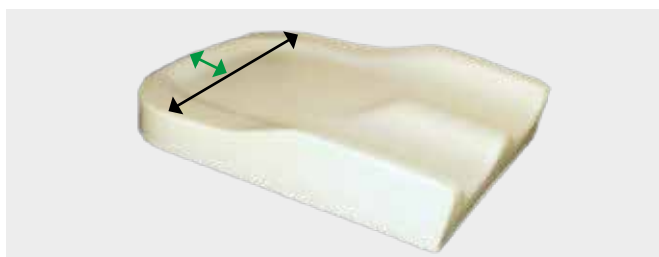

Rear cushion radius: +3.8 cm extra depth is designed to fill the gap between the seat and backrest/upholstery. **1.25 cm decrease in total width** allows cushion to fit both stated size and seats that are 2.5 cm narrower.

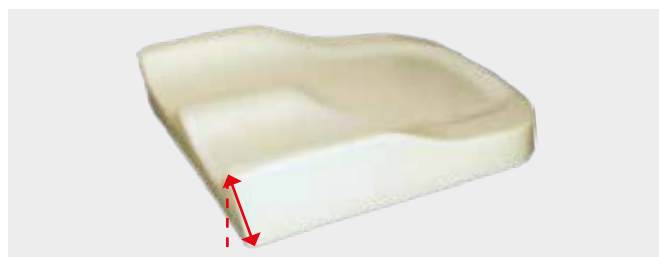

Waterfall front end

A diagonal cut along the front accommodates tight hamstrings, facilitates foot propulsion and standing transfers.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk


	Width	Depth	Height at adductor/abductor	Height at leg trough	Total product weight	Max. user weight	Warranty
KA-Vi	31 - 36 cm	31 - 36 cm	7,3 / 7 cm	5 cm	0,7 kg and up	68 kg (11 stone)	2 years
MA(WWDD)-Vi (Standard)	36, 41, 43, 46, 51, 53, 56 cm	38 - 51 cm	9,5 / 8,25 cm	6,35 cm	0,9 kg and up	136 kg (21 stone)	2 years

Colours and options

Black Meshtex

Kid*ab*ra fabric (Standard)

Pelvic Obliquity Kit

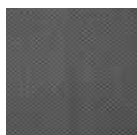

Easily installed on all Invacare *Matrx Vi Kid*ab*ra-Vi* cushions. Accommodates up to 25 mm obliquities